

sirius-XL ECLIPSE

TIMESCO™

Sirius-XL Eclipse blades feature a flexible tip designed for difficult intubations and are available in Macintosh 4, 3 and 2. Supplied as a blade only in a case or as a blade and handle set in a case.



CE marked



Fully autoclaveable



Latex Free



Extensive guarantees

**ISO
7376**

ISO 7376 compatible



Supplied with multilingual instructions



Compatible handles are colour coded

SPECIFICATION

FEATURES

Fully autoclaveable metal blade, covered by a 5 year guarantee or 4,000 autoclaves, whichever is sooner. See website for full terms and conditions.

Atraumatic distal end.

Tilting tip for elevation of the epiglottis, providing an improved view of the larynx, without excessive force.

Clear product information etched on the blade.


Green ISO 7376 fitting metal block compatible with fibre optic handles.

3 metal ball bearings for secure fitting.

Lever positioned at the rear of the blade to control the tilting tip.

The non-replaceable glass fibre bundle is fully enclosed in the blade which allows easy cleaning.

Rigid metal construction from fully autoclaveable, antimagnetic stainless steel.

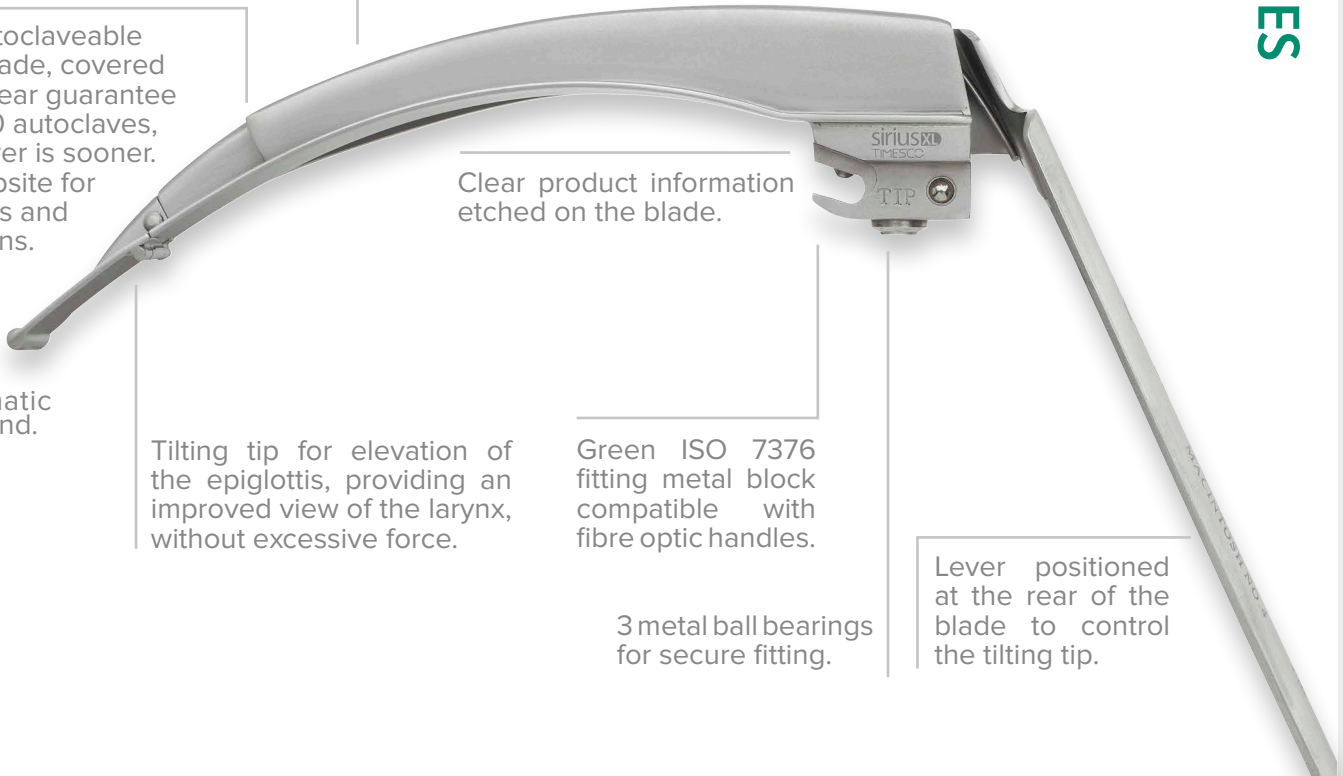



Product Code	Description	Measurements*
2955.175.15	Sirius XL Eclipse Macintosh Size 2	115mm x 24mm
2955.175.20	Sirius XL Eclipse Macintosh Size 3	135mm x 24mm
2955.175.25	Sirius XL Eclipse Macintosh Size 4	160mm x 24mm

Timesco Healthcare Ltd Timesco House 3 Carnival Park
Carnival Close Basildon Essex SS14 3WN England

T +44 (0) 1268 297 700 | F +44 (0) 1268 297 800
E info@timesco.com | timesco.com

L.L.058/1 E & O.E Timesco Feb 2018.

*Length x Height. Height is measured from the top of the block to the top of the blade. Measurements rounded to the nearest millimeter.