



XLED LARYNGOSCOPE HANDLES

NEW • IMPROVED • LONG LASTING • 3 x BRIGHTER • FLICKER FREE LED LAMP

TIMESCO
INSTRUMENTS FOR LIFE

New & improved with our best LED lamp yet...

RELIABLE ILLUMINATION...

Timesco's new and improved XLED laryngoscope handle is designed to outperform others and provide consistent, bright, reliable illumination when you need it most.

BRIGHT LIGHT, EVERY TIME...


Improved specification gives a bright, controlled illumination at all times. The lux and colour temperature from the lamp is consistent throughout the lifespan of the lamp and the batteries.

FLICKER FREE...

High performance and reliability from a clear, flicker free replaceable LED lamp. Intelligent electronics control the light intensity to be flicker free, even at low battery voltages.

NEW LAMP AT A GLANCE...

- Super-bright replaceable 3V LED lamp
- **Flicker free, white illumination**
- An exceptional lux output of 40,000 (at lamp)
- **Long lasting and improved bulb life**
- Colour temperature of 5,000(k)
- **3 x more lux output than other leading brands**

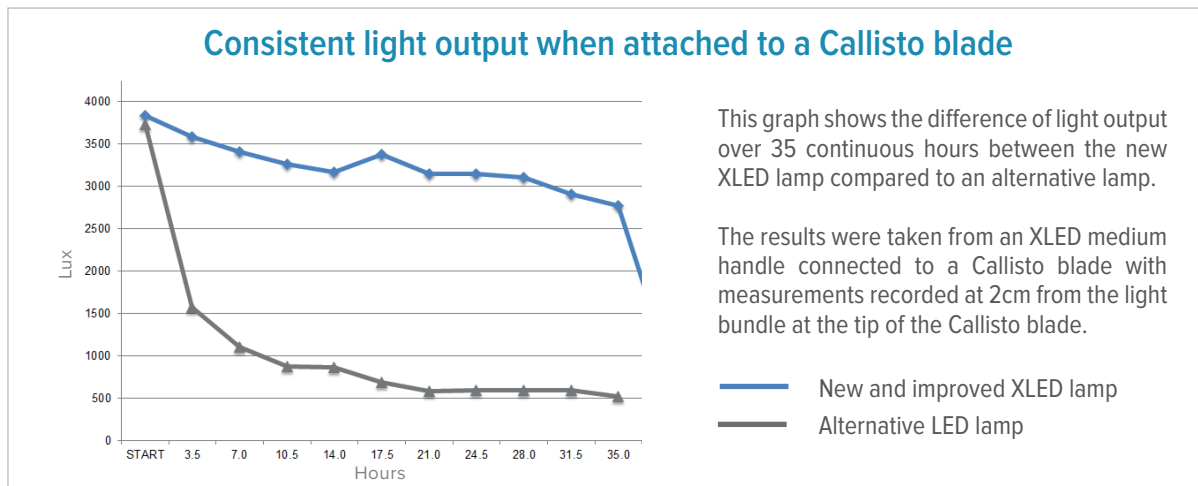

PRODUCT SPECIFICATION



Increased bulb lifespan



Latex free



Supplied with multilingual instructions



Fully autoclavable (excluding lamp)



ISO 7376 compatible



Extensive guarantee



Compatible blades are colour coded



Replaceable LED lamp



CE marked

PRODUCT FEATURES

The Timesco XLED handles provide exceptionally bright white illumination from a replaceable LED lamp supplying consistent, flicker free light & maximising battery life due to low power consumption. They are available in medium, small & stubby size and are compatible with reusable & single use blades.

A polished metal collar reflects and directs light up towards the fibre core of the blade for maximum illumination.

Super-bright replaceable 3V LED lamp providing excellent, flicker free, white illumination with a lux output of 40,000 (at lamp) and a colour temperature of 5,000(k).

Clear product branding with the product code and CE located on the opposite side to the XLED logo. The LOT number is shown on the top of the handle.

Colour coded band to clearly identify it is ISO 7376 fitting compatible with green fitting reusable and single use blades.

Knurled handle design provides excellent grip.

The new XLED lamp provides a brighter, whiter light and works with dry cell or rechargeable batteries without requiring a voltage of 3.5v.

Fully autoclavable metal handle (excluding lamp), covered by a 5 year guarantee or 4,000 autoclaves, whichever is sooner. See website for full terms and conditions.

Secure screw cap for easy access to batteries.



PRODUCT RANGE

Product Code	Description	Measurements*
3000.500.05	XLED Laryngoscope Handle, Small	146mm x 19mm
3000.500.10	XLED Laryngoscope Handle, Medium	146mm x 29mm
3000.508.10	XLED Laryngoscope Handle, Stubby	111mm x 32mm
3000.950.20	Spare XLED Lamps (Pack of 6)	-----

For more information please call +44 (0) 1268 297 710, email sales@timesco.com or visit www.timesco.com

PRODUCT FEATURES

A quick, simple and low cost solution to convert existing dry cell laryngoscope handles to a rechargeable version. Simply replace the dry cell batteries with the equivalent Timesco R engine.

Our Timesco R engines convert Timesco reusable handles to a rechargeable system. No special end caps/modifications required. (Excluding stubby handles).

The batteries use Lithium Ion technology to provide quick charging in an ION-R charger, (typically less than 2 hours for a full charge) and deliver consistent and reliable power.

Side plates can be removed to daisy chain 3 chargers together to receive power from a single power supply.



Only the engine is inserted into the charger, reducing the risk of cross contamination from a handle to the charger.

Timesco R engines can be fully charged in an ION-R charger in less than 2 hours.

The colour changing LED indicates the health of the battery and the charging cycle.

PRODUCT RANGE

Product Code	Description
3115.150.20	Timesco R Engine Small (AA Size), 2.5v
3115.150.50	Timesco R Engine Medium (C Size), 2.5v
3125.100.16	Timesco R Engine Charger Insert (AA Size)
3125.100.17	Timesco R Engine Charger Insert (C Size)
3125.100.15	ION-R Charger Only (excluding power supply)
3125.999.10	ION-R Single Charger Set - Including charger, power supply & both Timesco R inserts (AA and C)
3125.999.20	ION-R Dual Charger Set - Including 2 x chargers, 2 x power supplies, 2 x each Timesco R inserts (AA and C)

For more information please call +44 (0) 1268 297 710, email sales@timesco.com or visit www.timesco.com

Timesco Healthcare Ltd Timesco House 3 Carnival Park
Carnival Close Basildon Essex SS14 3WN England

T +44 (0) 1268 297 700 | F +44 (0) 1268 297 800
E info@timesco.com | timesco.com