

ET TUBES

TIMESCO™

An extensive range of Timesco Endotracheal Tubes with a high quality valve for reliable control of cuff inflation and pressure & a smooth soft tip reducing the potential of tracheal trauma during intubation.

A large variety of sizes and styles available to suit the needs of the patient, all at a competitive price.

SPECIFICATION



CE marked



Sterile (Ethylene Oxide)



Single Use



Latex Free

FEATURES




Murphy Eye:

Provides an alternate pathway for gas flow in case of occlusion as standard.

High quality valve:

For reliable control of cuff inflation and pressure.

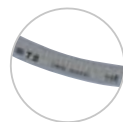

Closely spaced coil (on reinforced only):

Increases strength & flexibility, reducing the risk of occlusion whilst in use.

Smooth soft tip:

To reduce the potential of tracheal trauma during intubation.

High volume, low pressure cuff:

Minimising trauma to the trachea with a low pressure seal against the tracheal wall protecting the patient's airway, preventing leakage and allowing accurate control of gas flow.

Depth marker lines:

To facilitate the placement of the tube during intubation. (This should always be confirmed by auscultation).

Single use:

Reduces the risk of contamination.

Radio-opaque line:

To enable identification of the endotracheal tube on a chest X-Ray.

Calibration:

An accurate indication of insertion depth and an aid for monitoring the endotracheal tube position.



Product Type	Size Range Available	Box Quantity
Endotracheal Tube, Cuffed	3.0 to 10.0	20
Endotracheal Tube, Uncuffed	2.0 to 10.0	20
Endotracheal Tube, Preformed Nasal-North, Cuffed	3.0 to 8.5	20
Endotracheal Tube, Preformed Nasal-North, Uncuffed	3.0 to 8.5	20
Endotracheal Tube, Preformed Oral-South, Cuffed	3.0 to 8.5	20
Endotracheal Tube, Preformed Oral-South, Uncuffed	3.0 to 8.5	20
Endotracheal Tube, Reinforced, Cuffed	3.0 to 10.0	10
Endotracheal Tube, Reinforced, Uncuffed	3.0 to 10.0	10

For full product information visit www.timesco.com/et-tubes.

Please call 01268 297 710 or visit www.timesco.com for a quotation & further details.



Quality at a competitive price.

Our endotracheal tubes are available to buy now for up to 70% less than the average selling price of similar products on NHS Supply Chain helping you to save money whilst not compromising on quality.



End user satisfaction.

It's important that the end user is satisfied with the products we supply. We have ensured that our endotracheal tubes meet the standards & expectations for the performance of short duration tubes. Request a free sample by calling us.



Reputable brand.

As a supplier to the NHS and private sector with over 50 years' experience, Timesco are a reputable brand able to provide the best products and solutions for you.



Convenient logistics.

All of our products are available to order direct with our dedicated sales department who can be contacted via live chat, phone, email and fax.

Timesco Healthcare Ltd Timesco House 3 Carnival Park
Carnival Close Basildon Essex SS14 3WN England

T +44 (0) 1268 297 700 | F +44 (0) 1268 297 800

E info@timesco.com | timesco.com

L.L.065/1 E & O.E Timesco January 2019.

Image shown is Endotracheal Tube, Cuffed, Size 7 - TDM-D-6170