

**callisto** ECLIPSE  
TIMESCO™

The Timesco Callisto Eclipse blades are single use to reduce the risk of cross infection; with a rigid metal spatula, lever and tilting tip allowing greater flexibility and control during difficult intubations.

Available in Macintosh sizes 3 and 4 as blade only or as a convenient set with the Timesco Callisto Preloaded Handle.



CE marked

**ISO  
7376**

ISO 7376  
compatible



Single Use



Supplied with  
multilingual  
instructions



Latex Free



Compatible  
handles are  
colour coded


## SPECIFICATION

## FEATURES

Highly polished 4mm acrylic core for maximum illumination.

Supplied in boxes of 10 when sold individually or boxes of 5 when sold as a set.

Rigid metal spatula looks and feels like a reusable fibre optic blade.

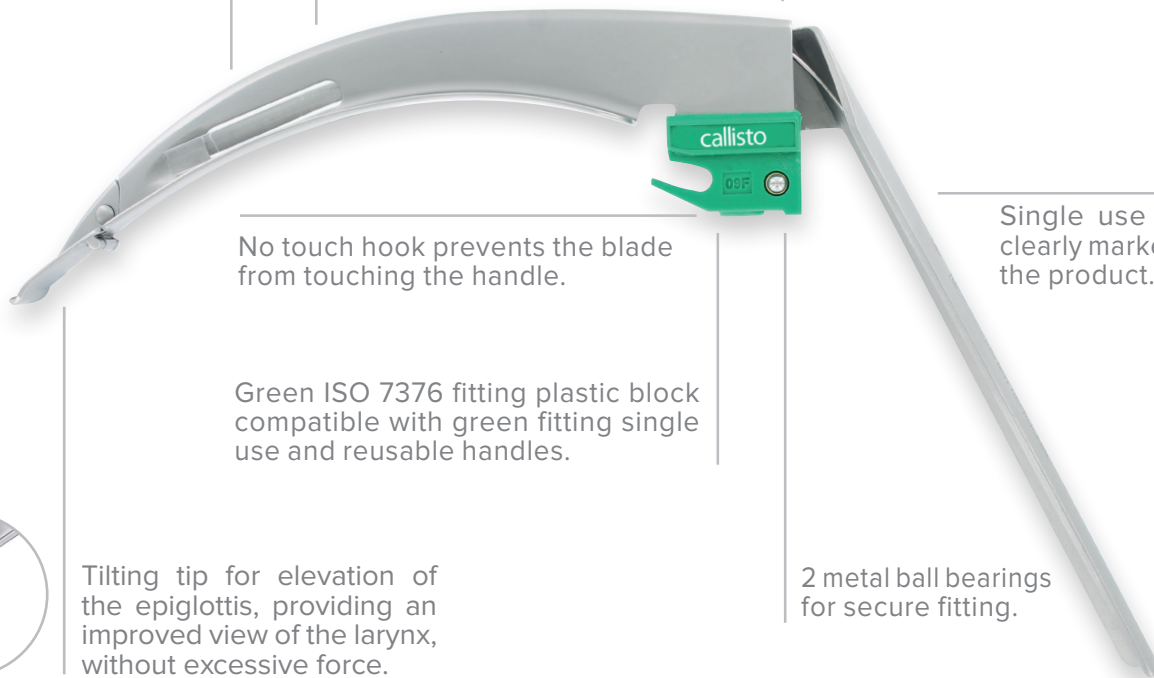

No touch hook prevents the blade from touching the handle.


Single use logo clearly marked on the product.

Green ISO 7376 fitting plastic block compatible with green fitting single use and reusable handles.

2 metal ball bearings for secure fitting.



Tilting tip for elevation of the epiglottis, providing an improved view of the larynx, without excessive force.

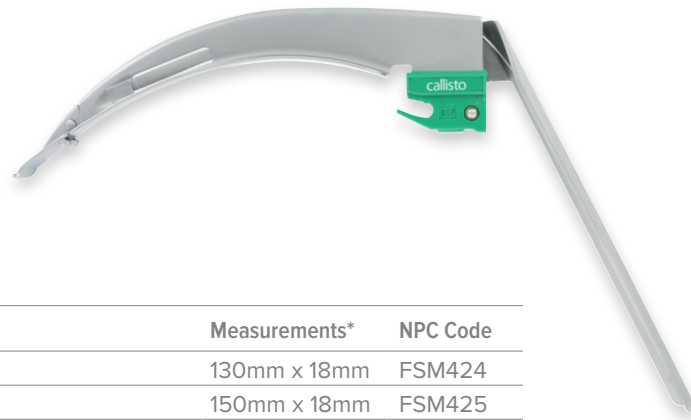

## CALLISTO ECLIPSE BLADES

Supplied in boxes of 10.

Product Code	Description	Measurements*	NPC Code
DS.3940.175.20	Callisto Eclipse Blade Size 3	130mm x 18mm	FSM424
DS.3940.175.25	Callisto Eclipse Blade Size 4	150mm x 18mm	FSM425

## CALLISTO ECLIPSE SETS

Supplied in boxes of 5.

Product Code	Description	Measurements of blade <sup>1</sup>	Measurements of handle <sup>2</sup>	NPC Code
DS.3950.953.20	Callisto Eclipse Blade Size 3 with Callisto Preloaded Handle (Medium)	130mm x 18mm	139mm x 33mm	FSM1835
DS.3950.953.25	Callisto Eclipse Blade Size 4 with Callisto Preloaded Handle (Medium)	150mm x 18mm	139mm x 33mm	FSM1836



Timesco Healthcare Ltd Timesco House 3 Carnival Park  
Carnival Close Basildon Essex SS14 3WN England

T +44 (0) 1268 297 700 | F +44 (0) 1268 297 800  
E [info@timesco.com](mailto:info@timesco.com) | [timesco.com](http://timesco.com)