



TIMESCO ENDOTRACHEAL TUBES

Timesco Endotracheal Tubes are available in cuffed and uncuffed options, featuring a high-quality valve for reliable cuff inflation control. Designed with a smooth, soft tip to reduce the risk of tracheal trauma.

 CE Marked

 Latex Free

 Single Use

 Sterilised using ethylene oxide


SPECIFICATION

FEATURES




Murphy Eye:
Provides an alternate pathway for gas flow in case of occlusion as standard.

High quality valve:
For reliable control of cuff inflation and pressure.

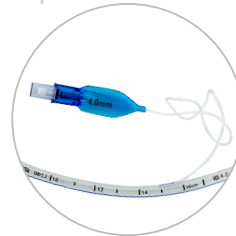

Radio-opaque line:
To enable identification of the endotracheal tube on a chest X-Ray.

Smooth soft tip:
To reduce the potential of tracheal trauma during intubation.

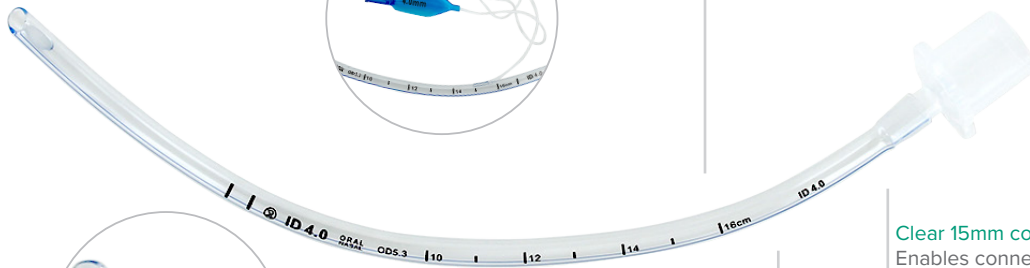

High volume, low pressure cuff:
Minimising trauma to the trachea with a low pressure seal against the tracheal wall protecting the patient's airway, preventing leakage and allowing accurate control of gas flow.




Single use:
Reduces the risk of contamination.

Depth marker lines:
Marker lines go all around the tube to facilitate its placement during intubation. (This should always be confirmed by auscultation).

Clear 15mm connector:
Enables connection to standard anaesthetic and ventilation breathing system equipment.

Calibration:
An accurate indication of insertion depth and an aid for monitoring the endotracheal tube position.

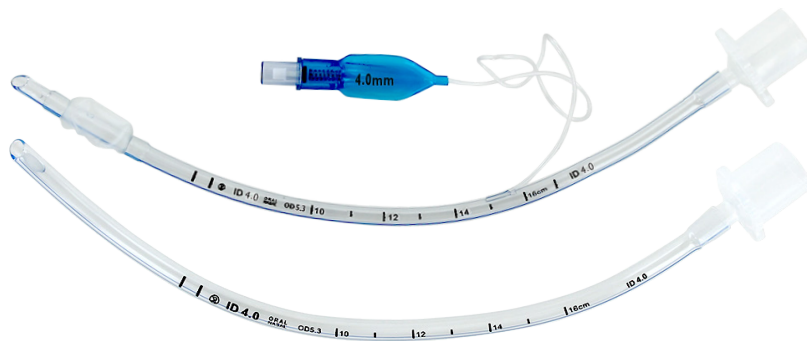

CAT. NO.	DESCRIPTION	NPC CODE
TDM-D-5130	Endotracheal Tube Cuffed 3.0 (Box of 20)	N/A
TDM-D-5135	Endotracheal Tube Cuffed 3.5 (Box of 20)	N/A
TDM-D-5140	Endotracheal Tube Cuffed 4.0 (Box of 20)	N/A
TDM-D-5145	Endotracheal Tube Cuffed 4.5 (Box of 20)	N/A
TDM-D-5150	Endotracheal Tube Cuffed 5.0 (Box of 20)	N/A
TDM-D-5155	Endotracheal Tube Cuffed 5.5 (Box of 20)	N/A
TDM-D-5160	Endotracheal Tube Cuffed 6.0 (Box of 20)	N/A
TDM-D-5165	Endotracheal Tube Cuffed 6.5 (Box of 20)	N/A
TDM-D-5170	Endotracheal Tube Cuffed 7.0 (Box of 20)	N/A
TDM-D-5175	Endotracheal Tube Cuffed 7.5 (Box of 20)	N/A
TDM-D-5180	Endotracheal Tube Cuffed 8.0 (Box of 20)	N/A
TDM-D-5185	Endotracheal Tube Cuffed 8.5 (Box of 20)	N/A
TDM-D-5190	Endotracheal Tube Cuffed 9.0 (Box of 20)	N/A
TDM-D-5195	Endotracheal Tube Cuffed 9.5 (Box of 20)	N/A
TDM-D-5199	Endotracheal Tube Cuffed 10.0 (Box of 20)	N/A
TDM-D-5220	Endotracheal Tube Uncuffed 2.0 (Box of 20)	N/A
TDM-D-5225	Endotracheal Tube Uncuffed 2.5 (Box of 20)	N/A
TDM-D-5230	Endotracheal Tube Uncuffed 3.0 (Box of 20)	N/A
TDM-D-5235	Endotracheal Tube Uncuffed 3.5 (Box of 20)	N/A
TDM-D-5240	Endotracheal Tube Uncuffed 4.0 (Box of 20)	N/A
TDM-D-5245	Endotracheal Tube Uncuffed 4.5 (Box of 20)	N/A
TDM-D-5250	Endotracheal Tube Uncuffed 5.0 (Box of 20)	N/A

Timesco Healthcare Ltd., 3 Carnival Park, Carnival Close,
 Basildon, Essex, SS14 3WN, United Kingdom
 T: +44 (0)1268 297 700 E: info@timesco.com W: timesco.com