




# optima view VIDEO LARYNGOSCOPE

FOR USE WITH REUSABLE BLADES

**TIMESCO**  
INSTRUMENTS FOR LIFE

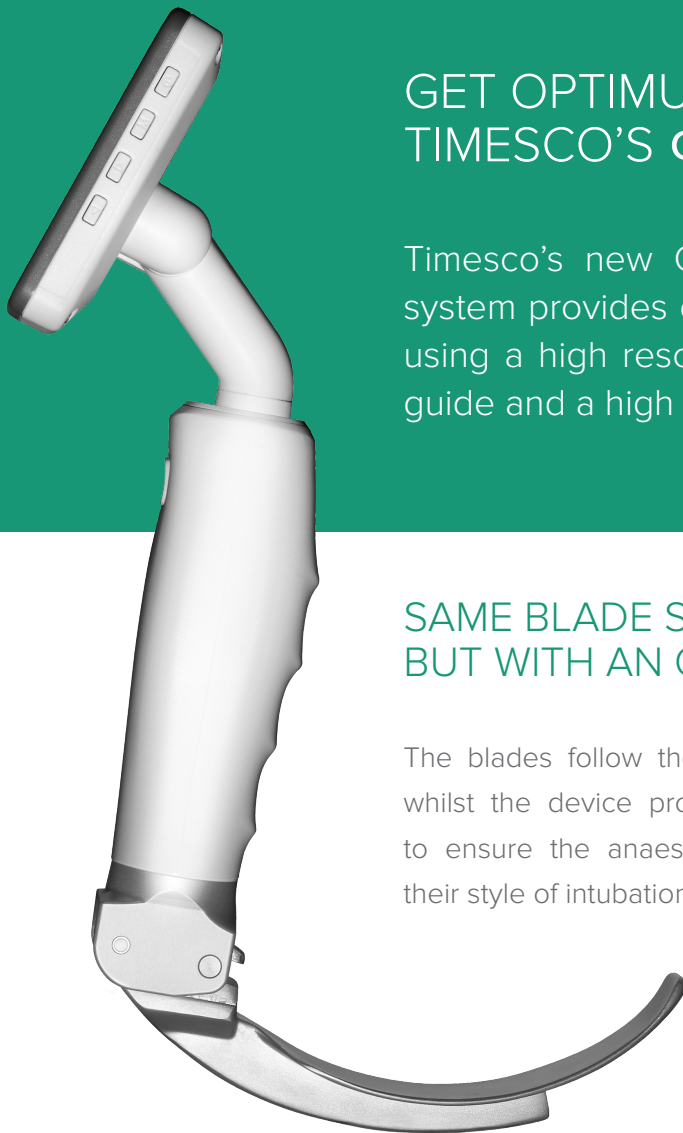

## GET OPTIMUM RESULTS WITH TIMESCO'S **optima view**

Timesco's new Optima View video laryngoscope system provides crystal clear colour correct outputs using a high resolution camera, a 5000k LED light guide and a high quality 3.75" LCD screen

## SAME BLADE SHAPE, SAME LIGHT OUTPUT BUT WITH AN OPTIMUM VIEW

The blades follow the same shape as a traditional Mac blade whilst the device provides a powerful white LED light output, to ensure the anaesthetist does not have to compromise on their style of intubation

### HIGH RESOLUTION

Large 3.75" LCD screen with a 2.0 megapixel blade camera

### AFFORDABLE

Optimum quality product at an affordable price point

### GUARANTEED

1 year guarantee and after-sales support available

Convenient rechargeable battery to allow ongoing use, lasting for over 2 hours

HDMI connectivity to external monitor or TV screen

Battery charge % indicator

Quick data transfer via USB

The screen rotates 180° from left to right and tilts 110° up and down. It also has a view angle of 160° making it convenient for teaching with multiple viewers

Blades are available in sizes Mac 1 - Mac 5 as well as Miller 00 - Miller 1. Shaped similarly to traditional Macintosh blades



The 2.0 megapixel, high resolution camera built into the blade, features a 50°-60° lens angle

Crystal clear colour correct image quality with a 3.75" LCD high resolution screen.  
Resolution ratio: 800 x 480RGB

Button Functions:  
1 On/Off  
2 Open menu/select  
3 Up scroll video/photo/menu  
4 Down scroll video/photo/menu/format disk (hold down)


24hr time and date display

Built-in 7.5GB storage, which can store approx. 10 hours of video or 70k+ pictures

One-click button suitably placed on the ergonomic handle for photography and video recording

Anti-fog camera lens built into the reusable blade

The reusable blade simply hooks and clips onto the Optima View system

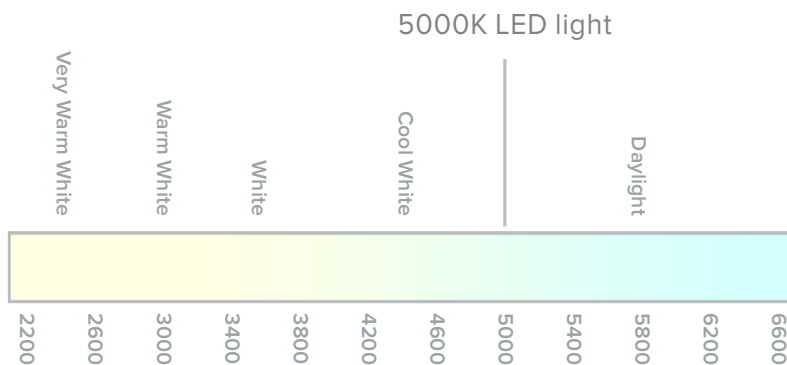



### Mac 1

Length: 101.5mm  
Width: 18.5mm  
Height: 34.4mm  
4500.150.10



### Mac 2

Length: 115mm  
Width: 21mm  
Height: 36.5mm  
4500.150.15



### Mac 3

Length: 126.8mm  
Width: 21.5mm  
Height: 39.5mm  
4500.150.20



### Mac 4

Length: 136.5mm  
Width: 21.5mm  
Height: 43mm  
4500.150.25



### Mac 5

Length: 149.7mm  
Width: 23mm  
Height: 50.8mm  
4500.150.30



### Miller 00

Length: 72.2mm  
Width: 14.6mm  
Height: 29.5mm  
4500.185.03



### Miller 0

Length: 84.3mm  
Width: 15mm  
Height: 29.5mm  
4500.185.05



### Miller 1

Length: 110mm  
Width: 16.7mm  
Height: 29.5mm  
4500.150.10



## PRODUCT VIDEO

Scan to watch product video.  
Alternatively, visit [timesco.com/optima-view/](https://timesco.com/optima-view/)



## FREE DEMO

Scan to organise a free demo.  
Alternatively, visit [timesco.com/optima-view/#demo](https://timesco.com/optima-view/#demo)

## PRODUCT SPECIFICATION

### Display

3.75" screen size

LCD monitor

800 x 480 RGB resolution ratio

4:3 aspect ratio

30 FPS video refresh rate

LED illumination

160° view angle

110° tilt angle

180° rotation angle

### Data output

Create & store files

### Internal memory

7.5GB internal storage

### Image output

USB interface to PC

Micro HDMI output

### Laryngoscope blade camera

2.0 megapixels resolution ratio

>800LUX illuminance

Video/image capture button on handle

### Includes

USB power adapter

USB charging cable

HDMI Cable

Hard case

### Handle

Video/image capture function

Comfortable ergonomic short handle design

### Battery

Rechargeable lithium battery

≥120 minutes on a full charge

3.7 volts

3200 mAh capacity

4-6 hours charging time

### Guarantee

Guarantee - 1 year for video device and blades

## PRODUCT RANGE

### Optima View video laryngoscope handle for use with reusable blades

4000.300.50	Optima View video laryngoscope handle and LCD display 3.75" for reusable blades
-------------	---

### Reusable Mac blades for use with Optima View

4500.150.10	Optima View Video Reusable Lary Mac 1 Blade
-------------	---

4500.150.15	Optima View Video Reusable Lary Mac 2 Blade
-------------	---

4500.150.20	Optima View Video Reusable Lary Mac 3 Blade
-------------	---

4500.150.25	Optima View Video Reusable Lary Mac 4 Blade
-------------	---

4500.150.30	Optima View Video Reusable Lary Mac 5 Blade
-------------	---

### Reusable Miller blades for use with Optima View

4500.185.03	Optima View Video Reusable Lary Miller 00 Blade
-------------	---

4500.185.05	Optima View Video Reusable Lary Miller 0 Blade
-------------	--

4500.185.10	Optima View Video Reusable Lary Miller 1 Blade
-------------	--

Classification: TIMESCO products are classified under the Harmonised Tariff: 90-1890-8400

## SINGLE-USE OPTIMA VIEW VIDEO LARYNGOSCOPE AND BLADES ALSO AVAILABLE

For more information please feel free to contact us:

T +44 (0) 1268 297 710

E [info@timesco.com](mailto:info@timesco.com)