



optima view VIDEO LARYNGOSCOPE

FOR USE WITH SINGLE-USE BLADES

TIMESCO
INSTRUMENTS FOR LIFE



GET OPTIMUM RESULTS WITH TIMESCO'S **optima view**

Timesco's new Optima View video laryngoscope system provides crystal clear colour correct outputs using a high resolution camera, a 5000k LED light guide and a high quality 3.75" LCD screen

SAME BLADE SHAPE, SAME LIGHT OUTPUT BUT WITH AN OPTIMUM VIEW

The blades follow the same shape as a traditional Mac blade whilst the device provides a powerful white LED light output, to ensure the anaesthetist does not have to compromise on their style of intubation

HIGH RESOLUTION

Large 3.75" LCD screen with a 2.0 megapixel camera

AFFORDABLE

Optimum quality product at an affordable price point

GUARANTEED

1 year guarantee and after-sales support available

Convenient rechargeable battery to allow ongoing use, lasting for over 2 hours


HDMI connectivity to external monitor or TV screen

Battery charge % indicator

Quick data transfer via USB

The screen rotates 180° from left to right and tilts 110° up and down. It also has a view angle of 160° making it convenient for teaching with multiple viewers

Blades are available in sizes Mac 1 - Mac 5 including channelled Mac 3 - Mac 5. Shaped similarly to traditional Macintosh blades

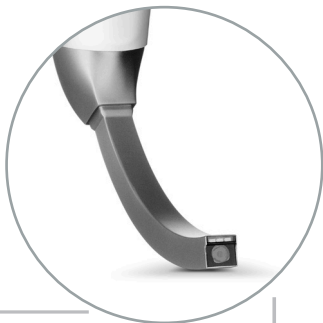

The 2.0 megapixel, high resolution camera features a 50°- 60° lens angle

Crystal clear colour correct image quality with a 3.75" LCD high resolution screen.
Resolution ratio: 800 x 480RGB

Button Functions:
1 On/Off
2 Open menu/select
3 Up scroll video/photo/menu
4 Down scroll video/photo/menu/format disk (hold down)

24hr time and date display


Built-in 7.5GB storage, which can store approx. 10 hours of video or 70k+ pictures

One-click button suitably placed on the ergonomic handle for photography and video recording

Anti-fog camera lens

The sterile, single-use blade simply slides and clips onto the Optima View system

5000K LED light

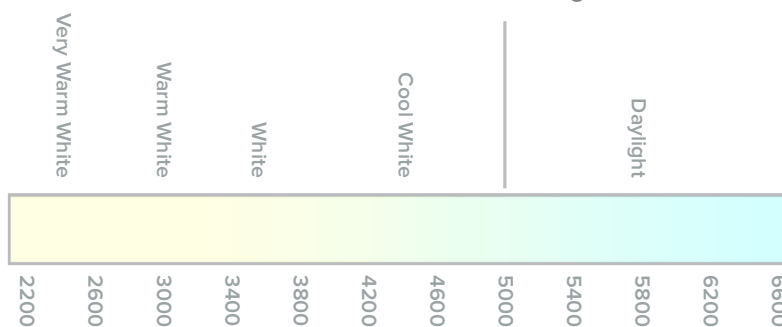



Mac 1

Length: 94mm
Width: 17.6mm
Height: 13mm

DS.4500.150.10



Mac 2

Length: 102mm
Width: 22mm
Height: 13mm

DS.4500.150.15



Mac 3

Length: 114mm
Width: 27mm
Height: 13mm

DS.4500.150.20



Channelled blade option available



Mac 4

Length: 127mm
Width: 29mm
Height: 13mm

DS.4500.150.25



Channelled blade option available



Mac 5

Length: 112mm
Width: 29mm
Height: 13mm

DS.4500.150.30



Channelled blade option available



PRODUCT VIDEO

Scan to watch product video.
Alternatively, visit timesco.com/optima-view/



FREE DEMO

Scan to organise a free demo.
Alternatively, visit timesco.com/optima-view/#demo

PRODUCT SPECIFICATION

Display

3.75" screen size

LCD monitor

800 x 480 RGB resolution ratio

4:3 aspect ratio

30 FPS video refresh rate

LED illumination

160° view angle

110° tilt angle

180° rotation angle

Data output

Create & store files

Internal memory

7.5GB internal storage

Image output

USB interface to PC

Micro HDMI output

Laryngoscope blade camera

2.0 megapixels resolution ratio

>800LUX illuminance

Video/image capture button on handle

Includes

USB power adapter

USB charging cable

HDMI Cable

Hard case

Handle

Video/image capture function

Comfortable ergonomic short handle design

Battery

Rechargeable lithium battery

≥120 minutes on a full charge

3.7 volts

3200 mAh capacity

4-6 hours charging time

Guarantee

Guarantee - 1 year for video device/handle

PRODUCT RANGE

Optima View video laryngoscope handle for use with single-use blades

4000.310.50	Optima View video laryngoscope handle and LCD display 3.75" for single-use blades
-------------	---

Single-use blades for use with Optima View

DS.4500.150.10	Optima View video laryngoscope blade. Mac 1 single-use sterile (pack of 10)
DS.4500.150.15	Optima View video laryngoscope blade. Mac 2 single-use sterile (pack of 10)
DS.4500.150.20	Optima View video laryngoscope blade. Mac 3 single-use sterile (pack of 10)
DS.4500.150.25	Optima View video laryngoscope blade. Mac 4 single-use sterile (pack of 10)
DS.4500.150.30	Optima View video laryngoscope blade. Mac 5 single-use sterile (pack of 10)

Single-use channelled blades for use with Optima View

DS.4500.151.20	Optima View video laryngoscope blade. Mac 3 single-use with channel sterile (pack of 10)
DS.4500.151.25	Optima View video laryngoscope blade. Mac 4 single-use with channel sterile (pack of 10)
DS.4500.151.30	Optima View video laryngoscope blade. Mac 5 single-use with channel sterile (pack of 10)

Classification: TIMESCO products are classified under the Harmonised Tariff: 90-1890-8400

REUSABLE OPTIMA VIEW VIDEO LARYNGOSCOPE AND BLADES ALSO AVAILABLE

For more information please feel free to contact us:

T +44 (0) 1268 297 710

E info@timesco.com